




INFEKTOCIDE SPRAY BR 503

Rapid Spray Disinfectant for medical devices

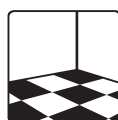

Fields of application

- Ideal for use in Hospital care and nursing homes, fitness and wellness studios, schools, nursery schools and health care facilities
- For all wet areas, sanitary facilities and for swimming pools as well as for sunbeds and fitness equipment
- For all water-resistant materials, surfaces and floors

Product advantages/properties

- Residue-free, thus no coating formed after application
- No wiping down or rinsing off necessary
- Acts quickly
- Ready to use
- Non-scented
- FDA Approved
- For disinfection of critical areas

Effective against Bacteria incl. MRSA & TBC, yeast, envelope viruses

1. Acinetobacter baumannii ATCC 1920
2. Bordetella pertusis ATCC
3. Capapenemase producer Klebsiella pneumoniae (hospital strain)
4. Escherichia coli ATCC 25922
5. ESBL Escherichia coli (Hospital strain)
6. Klebsiella pneumoniae ATCC 4352
7. Acinetobacter baumannii MDR strain (Hospital strain)
8. Meticillin Resistant Staphylococcus aureus (MRSA) (Hospital strain)
9. Meticillin Resistant Staphylococcus epidermidis (MRSE) (Hospital strain)
10. Proteus mirabilis ATCC 14153
11. Pseudomonas aeruginosa ATCC 10145
12. Salmonella typhi ATCC 6539
13. Staphylococcus aureus ATCC 6538
14. Streptococcus pneumoniae ATCC 49619
15. Streptococcus pyogenes (Hospital strain)
16. Vancomycin Resistant Staphylococcus aureus (VRE) (Hospital strain)
17. Clostridium difficile ATCC 43598
18. Candida albicans ATCC 10231

Technical data	
pH-Value	7.5 - 8.5
Product colour	Clear Color Liquid
Container*	5 L canister

* Please note: Some products may not be available in all territories. Please talk to your customer service representative.





INFEKTOCIDE SPRAY BR 503

Rapid Spray Disinfectant for medical devices

Application and dosage

- Spraying should be done with manual sprayer only.
- Spray till surface becomes visibly wet & allow it to dry for 10-15 min for optimum result.

For mopping / spraying:

- Switch off all electrical / electronic equipments before spraying them.
- Spray surface medical devices and lights etc. thoroughly from a distance of approx. 30cms.
- Spray 20-40 ml. of solutions per sq meter, for instant disinfection and allow alcohol to dry.
- Switch off all electrical / electronic equipment before spraying them.
- Spray surface medical devices and lights etc. thoroughly from a distance of approx. 30cms.
- Bactericides and levurocidal routine disinfection
: Pure - 1 min
- Noroviruses : Pure - 1min

Note

Use biocides safely. Always read the label and product information before use. Surfaces that come into contact with food or the skin must be rinsed with drinking water after disinfection. Absorb product which has been accidentally spilled with a liquid-binding material (sand, diatomaceous earth, universal binding agent). Do not let product enter the sewer system or bodies of water undiluted. Once packaging is completely empty, rinse out twice with water and recycle. A storage temperature of 5 to 25 °C

Hazard statements

H226 Flammable liquid and vapour.
H319 Causes serious eye irritation.
H336 May cause drowsiness or dizziness.

For further safety precautions, please refer to the safety data sheet for Infektocide Spray BR 503



The data in our product information sheets are based on our experience and careful investigations, and correspond to our present state of knowledge. We are, however, not able to examine and influence in detail the varying conditions for usage and materials, and those at premises where products are to be used. For this reason, we can only provide general information on usage. The quality of the work depends on expert assessment of the premises where products are to be used and how they are in fact employed by the user. In the case of doubt, the latter is obliged at all times to test the compatibility of the material to be treated with the Buzil Rossari product on an inconspicuous area or to obtain relevant advice on application technology. Liability for incomplete or incorrect data in our product information sheets shall only be assumed in cases of negligence caused by gross negligence or intent; any claims arising from the Product Liability Act shall remain unaffected. All previous product information shall lose its validity through the above listed product information.

Date 11.06.2019

BUZIL ROSSARI PRIVATE LIMITED

401, Omega, Hiranandani Business Park, Powai, Mumbai - 400076

T +91-2225701619 / 67411126 E brpl.admin@buzil-rossari.com W www.buzil-rossari.com